
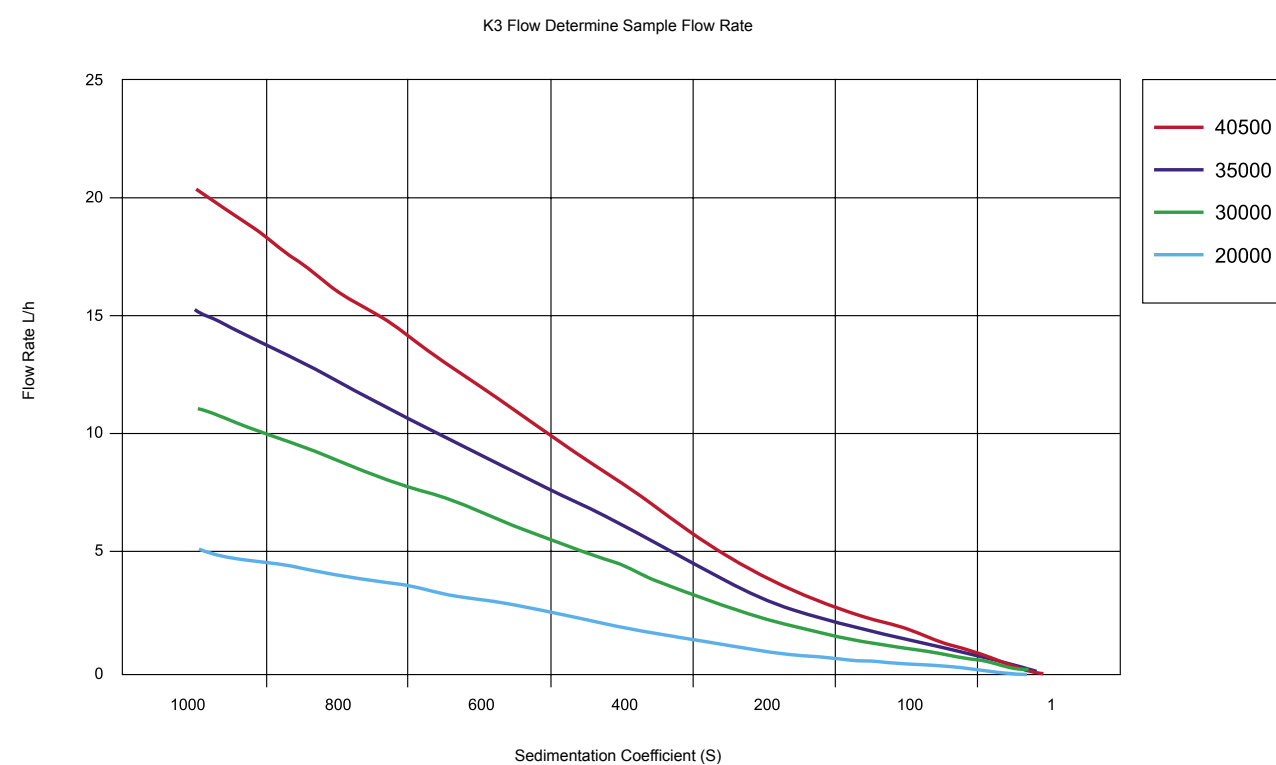


Alfa Wassermann

K3 Rotor for One-Step Solution, Concentration and Purification

Scaleable Rotors

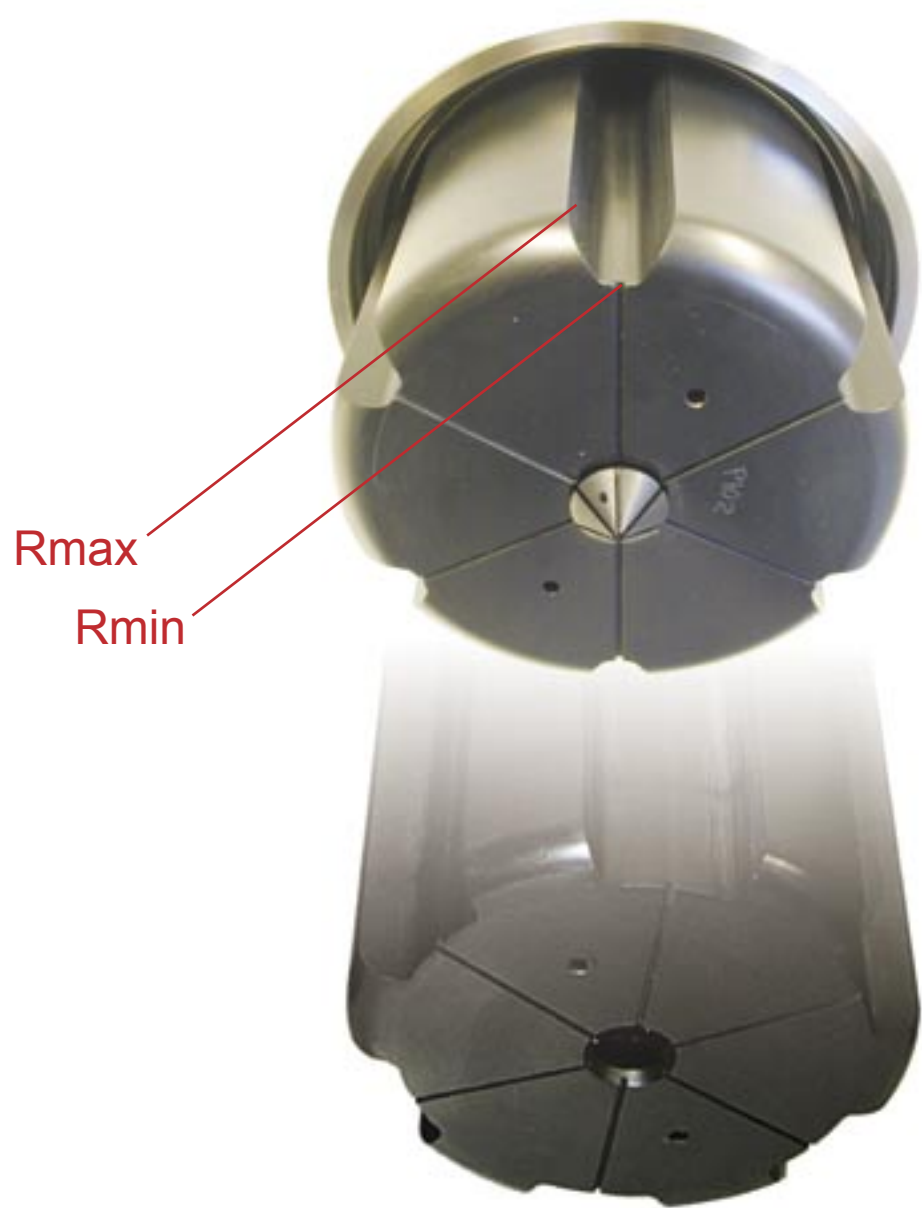

- Scaleable 16 fold
- 400ml to 3200ml rotors
- 40 500 rpm, 121 000 xg
- Up to 80 L/hr
- Titanium and Noryl® product contacts
- Low back-pressure flow channels

Typical Separation Profile

Speed: 40 500 rpm
 RCF: 121 200 xg
 Flow: 29 L/hr
 Gradient: 0-60% Sucrose
 Product Banding: 41% Sucrose
 Rotor: PK3
 Volume: 50 - 100 Litres
 Product Clean-out: 95%
 Product Recovery 70%
 Isopycnic banding process

